

Amazon 2025.01
Model GW0203



ULTRASONIC CLEANER

Instruction Manual >>

very noisy!

www.granbo.com

Professional teams provide design solutions for OEM/ODM

- Optional size, power and functions
- Frequency: 20/25/28/40/68/80/120/132/200 KHZ
- Voltage: Single-phase 110V/220V Three-phase 220V/380V/415V
- Material: SUS304/SUS316 T2.0/T2.5/T3.0
- Other functions: Filtration, Spray, Bubbling, Drying, Lifting, Throwing, Condensation, Slow-pulling, PLC (Programmable Logic Controller) system, etc.

Granbo Technology Industrial Shenzhen Co., Ltd.

Add: S&T Building, Quanzhi Technology Innovation Park, Shajing Street, Bao'an District, Shenzhen, China.

Tel: +86-755-27502790

E-mail: granbo@granbo.com

Website: www.granbo.com

www.granbosonic.com



APPLICATION:

Electronic factory, Car work shops, Industrial and mining area, Laboratories, Hospitals, Dental clinics, Watch shops, Optical shops, Jewelry shops, Mobile phone repair shops, and household use.

**Industrial Ultrasonic AIO machine**

- ① Capacity:14L-130L
- ② Extra-exciting Board, can keep work for 24hours
- ③ Power adjustable, with heating function

**Industrial Ultrasonic Single Tank machine**

- ① Capacity:40L-360L
- ② Extra-exciting Board, can keep work for 24hours
- ③ Power adjustable, with heating function

**Industrial Ultrasonic Double tank machine (Filtration cycle + drying)**

- ① Capacity:40L-360L
- ② Cleaning, rinsing, filtering, drying, spraying, lifting
- ③ Power adjustable, with heating function

**Industrial Ultrasonic Multi-tank machine (clean +spray +cycle)**

- ① Support Multi-tank customized
- ② Cleaning, rinsing, filtering, drying, spraying, lifting
- ③ Power adjustable, with heating function



Small three-tank cleaning machine

- ① Capacity:3X6.5L
- ② Digital degassing + heating
- ③ Cleaning+ Cleaning+ Spraying




Vinyl record ultrasonic cleaner + liftable bracket

- ① Capacity:6.5L
- ② Wash up to 8pcs vinyl records at a time
- ③ With lifting function



Submersible Ultrasonic Processor

- ① Length, shape, voltage and material support customization
- ② Cleaning, decomposition, defoaming, mixing, emulsification, homogenization, extraction, catalysis, stirring, etc.
- ③ Power adjustable

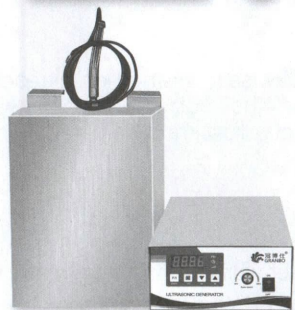

Submersible Ultrasonic Plate

- ① Support customization
- ② Bottom-mounted, side-mounted, top-mounted
- ③ Power adjustable

intro

Ultrasonic cleaner

- Thanks for purchasing this machine.
- Read all instructions before operating the device.
- Please keep the manual for reference.

Catalog [Index]

• Security precautions & Daily maintenance & Application description	1-3	1
• Product introduction & Product specification & Operation steps	3-15	12
• Creative development & Non-standard customization	16-17	16

Model GW0203

Cautions

- Ensure the power supply is in rated range before power cable is connected. Refitting is strictly prohibited! Attention that the control panel will be corroded by organic solution, strong acid and strong alkali.
- Ensure the earth wiring is well connected before starting.
- Activate key of system should be on position of "ON/OFF" button.
- Do not turn on machine if the tank is empty or the transducer will be damaged.
- Please close the lid to reduce noise and pay attention to water and steam in case of burning skin while open the lid.
- Do not relocate the machine when fluid in the tank in case of overflow.
- Water-soluble liquid is suggested for bench top ultrasonic cleaner. Strong acid or flammability cleaner is strictly forbidden. If solvent type cleaner is added in tank, please stop working for a while after 30mins cycle.
- Do not use the machine in below environment.
 - The place where temperature change fiercely
 - The place where humidity is too high and is easy to produce dew.
 - The place where vibration or impact is strong
 - The place where exists corrosive gas or dust
 - The place where water, oil or chemicals splash
 - The place where is filled with explosive and flammable gas



Do not work for long time when you use flammable liquid like alcohol!
You will get better cleaning effect if cleaning solution is added in ultrasonic cleaner.

Maintenance

- The unit must be opened by authorised specialised person only for maintenance and care of it.
- Clean the contaminant in the tank frequently.

Trouble clearing

NO	Problems	Possible causes	Solutions	Remarks
1	No ultrasonic	A.Power supply not connect B.Fuse broken C.Cable short circuit D.Transducer short circuit E.PCB Board broken F.Other reasons	A.Check and plug power switch. B.Check the fitted power supply and change the same specification fuse. C.Connect fitted cable or replace a new. D.Inquiry our after service engineer. E.Check the broken part and replace it. F.Inquiry our after service engineer.	
2	Time control failure	A.Timer knob out of control B.Timer failure	A.Loosen or tighten the screw. B.Replace timer or digital panel. C.Inquiry our after service engineer.	
3	No heating	A.Heater power switch linkage. B.Fuse burn outl. C.Heating pad burn out. D.Digital display board out of control E.Other reasons	A.Check heating plug and connect well, Check outlet line with multimeter: if OK and resistance value is few hundred OHMs. B.Change the same specification fuse. C.Replace the bad heating pad if open circuit. D.Check the broken part and replace it. E.Inquiry our after service engineer.	Suggestion 50-60°C
4	Temperature control failure	A.Thermostat loosen. B.Thermostat tube broken. C.Digital display of control.	A.Fasten the thermostat header. B.Replace thermostat. C.Check the broken part and replace it. D.Inquiry our after service engineer.	
5	Not well cleaning	A.Not Strong ultrasonic cleaning. B.Too high too low liquid surfacel. C.too high too low temperature. D.Not suitable cleaning liquid. E.Other reasons	A.Connect ultrasonic button and adjust. B. Adjust liquid into best urface. C.Adjust temperature into the most fitted. D.Stop and switch off power supply, replace suitable liquid after the previous liquid cool down. E.Inquiry our after service engineer.	Suggestion 50-60°C
6	Electric leakage	C.Customer side not rounded. D.Machine not grounded.	C.Make sure the wire. D.Check if machine earth loosen.	
7	Other problems			

X Series (Timer & Heater & Degas & Semiwave & Adjustable power) operation steps



Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.



Uncover the lip and fill the tank with tap water.

No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water, but do not over the MAX line.



Time setting

Default time is "05 00". Press "TIME+" once to increase time by 1minute.Press "TIME-" once to decrease time by 1minute. Press ON/OFF after time setting and ultrasonic indicator light on, machine works follow the setting parameter and stop work when time counts down to "00 00", then indicator light off.



Temperature

Temperature display panel shows "50" and actual detergent temp. Press "TEMP+" to increase setting temperature by 1°C, press "TEMP-" to decrease setting temperature by 1°C. Press "ON/OFF" to start heating function and heater indicator on, actual temperature number will blazing by 1Hz Frequency during heating process, heater indicator light will off when temperature reaches the setting temperature.



Degas and Semiwave

After time setted , press "DEGAS" button to switch into degassing mode , then the machine will work alternately by semi-wave for 5s then full wave for 5s , it can exhaust air and prevent oxidation. When device on standby, press "SEMIWAVE" button, switch into semiwave working, the ultrasonic power output reduce by half, but the ultrasound strength only reduce by 20% , this mode good for fine cleaning and low power operation requirements.



Power


Turn the power adjustment knob clockwise to strengthen ultrasonic power, maximum is rated power. Turn the knob counterclockwise to lower power, minimum is 0W. Do not adjust knob frequently when device working.



Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.

X Series Product introduction

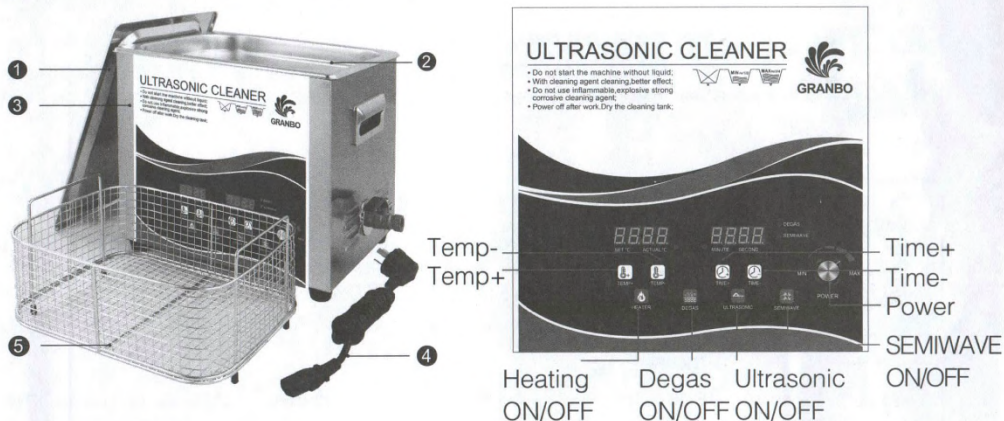

- ① Lid ③ Stainless steel body ⑤ Stainless steel basket
② Tank ④ Plug

Product specification

Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GX0203	3.2	40	0-120	150	240x135x100	265x165x230	20-80	0-30
GX0303	3.2	40	0-180	150	240x135x100	265x165x230	20-80	0-30
GX0304	4.5	40	0-180	300	300x150x100	325x175x230	20-80	0-30
GX0306	6.5	40	0-180	300	300x150x150	325x175x270	20-80	0-30
GX0410	10	40	0-240	300	300x240x150	325x265x275	20-80	0-30
GX0610	10	40	0-360	300	300x240x150	325x265x275	20-80	0-30
GX0615	15	40	0-360	450	330x300x150	355x325x275	20-80	0-30
GX0815	15	40	0-480	450	330x300x150	355x325x275	20-80	0-30
GX0915	15	40	0-540	450	330x300x150	355x325x275	20-80	0-30
GX0822	22	40	0-480	600	500x300x150	530x330x280	20-80	0-30
GX1022	22	40	0-600	600	500x300x150	530x330x280	20-80	0-30
GX1222	22	40	0-720	600	500x300x150	530x330x280	20-80	0-30
GX1422	22	40	0-840	600	500x300x150	530x330x280	20-80	0-30
GX1522	22	40	0-900	600	500x300x150	530x330x280	20-80	0-30
GX1030	30	40	0-600	600	500x300x200	530x330x330	20-80	0-30
GX1230	30	40	0-720	600	500x300x200	530x330x330	20-80	0-30
GX1430	30	40	0-840	600	500x300x200	530x330x330	20-80	0-30
GX1530	30	40	0-900	600	500x300x200	530x330x330	20-80	0-30


applications

Ultrasonic application description

- The standard ultrasonic cleaners adopt high quality stainless steel plate with strong corrosion resistance and long lifespan, transducers with advanced connection sticking technology, high electroacoustic conversion efficiency and strong ultrasonic output, also equipped with automatic constant temperature heating device. Temperature range: 20-80°C.
- Suit for varies industries and objects cleaning, such as precision electronic parts, eyeglasses & clock parts, optical glass parts, hardware, jewelry, semiconductor silicon wafer, printer inkjet, medical instruments and electroplate parts etc.

Industry	The cleaning products and materials	Clear dirt
Semi-conductor	Integrated circuit, power tube, silicon wafer, diode, lead frame, capillary, tray, etc.	Hards, etching oil, polishing wax, dust particles, etc.
Electrical & electronic machine	Tube parts, cathode ray tube, printed circuit board, quartz parts, electronic components, telephone switching equipment, speaker components, power meter, LCD glass, core iron parts, computer floppy disk, video parts, head, photo die mask, etc.	Finger print, powder, cutting oil, stamping oil, iron filings, polish-ing materials, walnut powder, polishing wax, resin, dust, etc.
Precision machine	Bearing, sewing machine parts, typewriter, textile, machine, optical mechanical device, gas valve, watches cameras, metal filter.	Machine cutting oil, iron filings, olishing powder, finger print, oil, grease, dirt, etc.
Optical device	Glasses, lens, prism, optical lens, filter lens, glass device, film, optical fiber, etc.	Plastic resin, paraffin, finger printing, etc.
Hardware & machinery parts	Bearing, gear, ball, metal shaft parts, tools, adjustable valve and cylinder parts, burner, compressors, hydraulic press, gun and ultracentrifuge, city water faucet, etc.	Cutting oil, iron filings, grease, polishing powder, finger printing and so on.
Medical instrument	Medical instrument, denture, etc.	Iron filings. Polishing powder, oil, stamping oil, dirt, etc.
Electroplate	Galvanized parts, mold, stamping parts, etc.	Polishing scrap iron, oil, black iron shell, rust, oxidation shell, scrap iron, polishing powder, stamping oil, dirt, etc.
Car parts	Piston ring, carburetor, flow meter housing, compressor shell, electrical components, etc.	Iron filings, polishing powder, oil stamping oil, dirt etc.
Chemical fiber	Chemical or artificial fiber nozzle filter protector chemical fiber texture, etc.	Chemical colloid, glue, and other solid material. dust. Etc.

B Series Product introduction

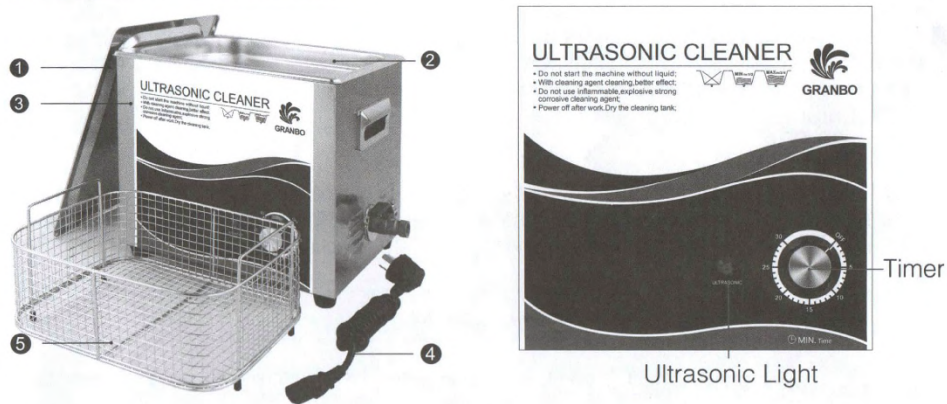

- ① Lid ③ Stainless steel body ⑤ Stainless steel basket
② Tank ④ Plug

Product specification

Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GB0101	1.3	40	60	-----	150x135x65	175x165x190	-----	0-30
GB0201	1.3	40	120	-----	150x135x65	175x165x190	-----	0-30
GB0102	2.0	40	60	-----	150x135x100	175x165x225	-----	0-30
GB0202	2.0	40	120	-----	150x135x100	175x165x225	-----	0-30
GB0203	3.2	40	120	-----	240x135x100	265x165x230	-----	0-30
GB0303	3.2	40	180	-----	240x135x100	265x165x230	-----	0-30
GB0304	4.5	40	180	-----	300x150x100	325x175x230	-----	0-30
GB0306	6.5	40	180	-----	300x150x150	325x175x270	-----	0-30
GB0410	10	40	240	-----	300x240x150	325x265x275	-----	0-30
GB0610	10	40	360	-----	300x240x150	325x265x275	-----	0-30
GB0615	15	40	360	-----	330x300x150	355x325x275	-----	0-30
GB0815	15	40	480	-----	330x300x150	355x325x275	-----	0-30
GB0915	15	40	540	-----	330x300x150	355x325x275	-----	0-30
GB0822	22	40	480	-----	500x300x150	530x330x280	-----	0-30
GB1022	22	40	600	-----	500x300x150	530x330x280	-----	0-30
GB1222	22	40	720	-----	500x300x150	530x330x280	-----	0-30
GB1422	22	40	840	-----	500x300x150	530x330x280	-----	0-30
GB1522	22	40	900	-----	500x300x150	530x330x280	-----	0-30
GB1030	30	40	600	-----	500x300x200	530x330x330	-----	0-30
GB1230	30	40	720	-----	500x300x200	530x330x330	-----	0-30
GB1430	30	40	840	-----	500x300x200	530x330x330	-----	0-30
GB1530	30	40	900	-----	500x300x200	530x330x330	-----	0-30

W Series (Timer & Heater & Degas & Semiwave) operation steps



Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.



Uncover the lip and fill the tank with tap water.

No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water, but do not over the MAX line.



Time setting

Default time is "05 00". Press "TIME+" once to increase time by 1minute. Press "TIME-" once to decrease time by 1minute. Press ON/OFF after time setting and ultrasonic indicator light on, machine works follow the setting parameter and stop work when time counts down to "00 00", then indicator light off.



Temperature

Temperature display panel shows "50" and actual detergent temp. Press "TEMP+" to increase setting temperature by 1°C, press "TEMP-" to decrease setting temperature by 1°C. Press "ON/OFF" to start heating function and heater indicator on, actual temperature number will blazng by 1Hz Frequency during heating process, heater indicator light will off when temperature reaches the setting temperature.



Degas

After time settled, press "DEGAS" button to switch into degassing mode, then the machine will work alternately by semi-wave for 5s then full wave for 5s, it can exhaust air and prevent oxidation.



Semiwave

After time settled, press "ULTRASONIC" button to switch into full wave working. Press "SEMIWAVE" button, device switch into half wave working, ultrasonic power output reduce by half, but the ultrasound strength only reduce by 20%, this mode good for fine cleaning and low power operation requirements. (1.3L and 2L models, press "DEGAS" button 3times within 3s to turn on semiwave mode.)



Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.

W Series Product introduction



- 1 Lid 3 Stainless steel body 5 Stainless steel basket
2 Tank 4 Plug

Product specification

Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GW0101	1.3	40	30/60	150	150x135x65	175x165x190	20-80	0-30
GW0201	1.3	40	60/120	150	150x135x65	175x165x190	20-80	0-30
GW0102	2.0	40	30/60	150	150x135x100	175x165x225	20-80	0-30
GW0202	2.0	40	60/120	150	150x135x100	175x165x225	20-80	0-30
GW0203	3.2	40	60/120	150	240x135x100	265x165x230	20-80	0-30
GW0303	3.2	40	90/180	150	240x135x100	265x165x230	20-80	0-30
GW0304	4.5	40	90/180	300	300x150x100	325x175x230	20-80	0-30
GW0306	6.5	40	90/180	300	300x150x150	325x175x270	20-80	0-30
GW0410	10	40	120/240	300	300x240x150	325x265x275	20-80	0-30
GW0610	10	40	180/360	300	300x240x150	325x265x275	20-80	0-30
GW0615	15	40	180/360	450	330x300x150	355x325x275	20-80	0-30
GW0815	15	40	240/480	450	330x300x150	355x325x275	20-80	0-30
GW0915	15	40	270/540	450	330x300x150	355x325x275	20-80	0-30
GW0822	22	40	240/480	600	500x300x150	530x330x280	20-80	0-30
GW1022	22	40	300/600	600	500x300x150	530x330x280	20-80	0-30
GW1222	22	40	360/720	600	500x300x150	530x330x280	20-80	0-30
GW1422	22	40	420/840	600	500x300x150	530x330x280	20-80	0-30
GW1522	22	40	450/900	600	500x300x150	530x330x280	20-80	0-30
GW1030	30	40	300/600	600	500x300x200	530x330x330	20-80	0-30
GW1230	30	40	360/720	600	500x300x200	530x330x330	20-80	0-30
GW1430	30	40	420/840	600	500x300x200	530x330x330	20-80	0-30
GW1530	30	40	450/900	600	500x300x200	530x330x330	20-80	0-30

B Series Mechanical model (Without heater) operation steps



Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.



Uncover the lip and fill the tank with tap water.


No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water, but do not over the MAX line.



Time setting

Clockwise rotate the time switch to chose the time you required between 0-30min. When the indicator light is on and make "ZIZI" sound, it shows ultrasonic works well.


Please take attention that anticlockwise rotate is forbidden.

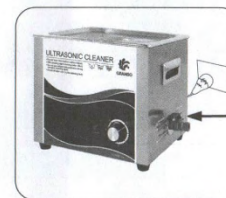

Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.

D Series Product introduction

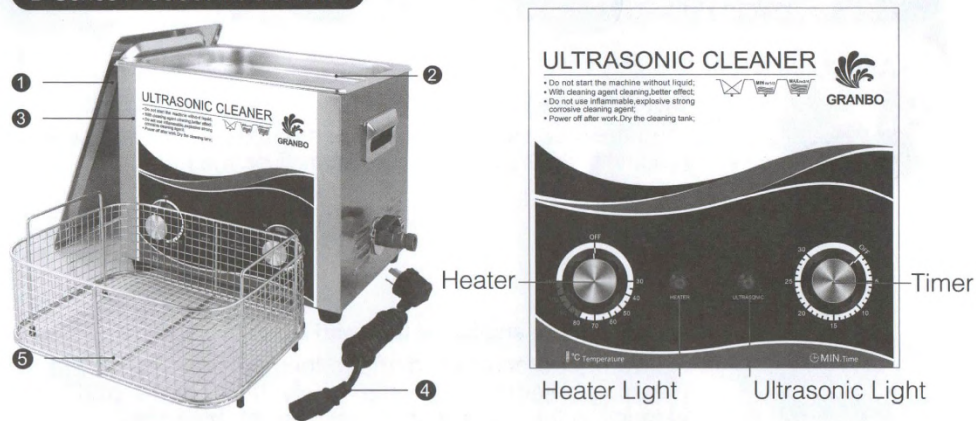

- 1 Lid 3 Stainless steel body 5 Stainless steel basket
2 Tank 4 Plug

Product specification

Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GD0101	1.3	40	60	150	150x135x65	175x165x190	20-80	0-30
GD0201	1.3	40	120	150	150x135x65	175x165x190	20-80	0-30
GD0102	2.0	40	60	150	150x135x100	175x165x225	20-80	0-30
GD0202	2.0	40	120	150	150x135x100	175x165x225	20-80	0-30
GD0203	3.2	40	120	150	240x135x100	265x165x230	20-80	0-30
GD0303	3.2	40	180	150	240x135x100	265x165x230	20-80	0-30
GD0304	4.5	40	180	300	300x150x100	325x175x230	20-80	0-30
GD0306	6.5	40	180	300	300x150x150	325x175x270	20-80	0-30
GD0410	10	40	240	300	300x240x150	325x265x275	20-80	0-30
GD0610	10	40	360	300	300x240x150	325x265x275	20-80	0-30
GD0615	15	40	360	450	330x300x150	355x325x275	20-80	0-30
GD0815	15	40	480	450	330x300x150	355x325x275	20-80	0-30
GD0915	15	40	540	450	330x300x150	355x325x275	20-80	0-30
GD0822	22	40	480	600	500x300x150	530x330x280	20-80	0-30
GD1022	22	40	600	600	500x300x150	530x330x280	20-80	0-30
GD1222	22	40	720	600	500x300x150	530x330x280	20-80	0-30
GD1422	22	40	840	600	500x300x150	530x330x280	20-80	0-30
GD1522	22	40	900	600	500x300x150	530x330x280	20-80	0-30
GD1030	30	40	600	600	500x300x200	530x330x330	20-80	0-30
GD1230	30	40	720	600	500x300x200	530x330x330	20-80	0-30
GD1430	30	40	840	600	500x300x200	530x330x330	20-80	0-30
GD1530	30	40	900	600	500x300x200	530x330x330	20-80	0-30

T Series (Timer & Heater & Degas & Adjustable power) operation steps



Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.



Uncover the lip and fill the tank with tap water.

No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water, but do not over the MAX line.



Time setting

Default time is "05 00". Press "TIME+" once to increase time by 1minute. Press "TIME-" once to decrease time by 1minute. Press ON/OFF after time setting and ultrasonic indicator light on, machine works follow the setting parameter and stop work when time counts down to "00 00", then indicator light off.



Temperature

Temperature display panel shows "50" and actual detergent temp. Press "TEMP+" to increase setting temperature by 1°C, press "TEMP-" to decrease setting temperature by 1°C. Press "ON/OFF" to start heating function and heater indicator on, actual temperature number will blazing by 1Hz Frequency during heating process, heater indicator light will off when temperature reaches the setting temperature.



Degas

In the downtime, press "ULTRASONIC" button 3 times within 3 seconds to enter the intermittent degassing mode, then the machine will work for 6 seconds and stop for 3 seconds to run alternately until the time end.



Power


Turn the power adjustment knob clockwise to strengthen ultrasonic power, maximum is rated power. Turn the knob counterclockwise to lower power, minimum is 0W. Do not adjust knob frequently when device working.



Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.

T Series Product introduction

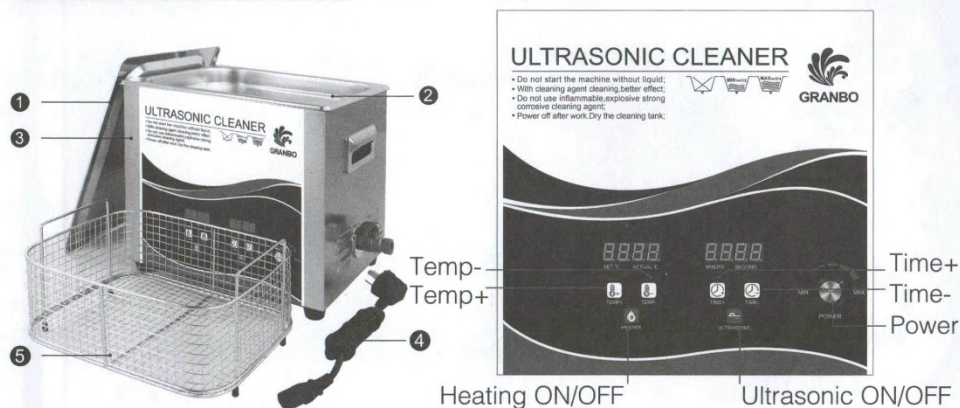

- 1 Lid 2 Tank 3 Stainless steel body 4 Plug 5 Stainless steel basket

Product specification

Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GT0101	1.3	40	0-60	150	150x135x65	175x165x190	20-80	0-30
GT0201	1.3	40	0-120	150	150x135x65	175x165x190	20-80	0-30
GT0102	2.0	40	0-60	150	150x135x100	175x165x225	20-80	0-30
GT0202	2.0	40	0-120	150	150x135x100	175x165x225	20-80	0-30
GT0203	3.2	40	0-120	150	240x135x100	265x165x230	20-80	0-30
GT0303	3.2	40	0-180	150	240x135x100	265x165x230	20-80	0-30
GT0304	4.5	40	0-180	300	300x150x100	325x175x230	20-80	0-30
GT0306	6.5	40	0-180	300	300x150x150	325x175x270	20-80	0-30
GT0410	10	40	0-240	300	300x240x150	325x265x275	20-80	0-30
GT0610	10	40	0-360	300	300x240x150	325x265x275	20-80	0-30
GT0615	15	40	0-360	450	330x300x150	355x325x275	20-80	0-30
GT0815	15	40	0-480	450	330x300x150	355x325x275	20-80	0-30
GT0915	15	40	0-540	450	330x300x150	355x325x275	20-80	0-30
GT0822	22	40	0-480	600	500x300x150	530x330x280	20-80	0-30
GT1022	22	40	0-600	600	500x300x150	530x330x280	20-80	0-30
GT1222	22	40	0-720	600	500x300x150	530x330x280	20-80	0-30
GT1422	22	40	0-840	600	500x300x150	530x330x280	20-80	0-30
GT1522	22	40	0-900	600	500x300x150	530x330x280	20-80	0-30
GT1030	30	40	0-600	600	500x300x200	530x330x330	20-80	0-30
GT1230	30	40	0-720	600	500x300x200	530x330x330	20-80	0-30
GT1430	30	40	0-840	600	500x300x200	530x330x330	20-80	0-30
GT1530	30	40	0-900	600	500x300x200	530x330x330	20-80	0-30

D Series Mechanical model (Timer & Heater) operation steps




Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.



Uncover the lip and fill the tank with tap water.

No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water, but do not over the MAX line.

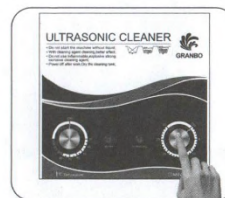

Time setting

Clockwise rotate the time switch to chose the time you required between 0-30min. When the indicator light is on and make "ZIZI" sound, it shows ultrasonic works well.

Please take attention that anticlockwise rotate is forbidden.



Heater setting

Rotate temperature control knob to adjusted temp and start heating, 40°C-60°C is suggested. (Heater is not suggesting for normal environment). When actual temp is lower than setting temp, heating function will auto start and the heat indicator light on. When temp reaches the setting temp, heat indicator light will off.



Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.

S Series Product introduction



- 1 Lid 3 Stainless steel body 5 Stainless steel basket
2 Tank 4 Plug

Product specification


Model	Capacity (L)	Frequency (KHz)	Ultrasonic Power (W)	Heating Power (W)	Tank Size (LxWxH mm)	Unit Size (LxWxH mm)	Temp. (°C)	Timer (Min)
GS0101	1.3	40	60	150	150x135x65	175x165x190	20-80	0-30
GS0201	1.3	40	120	150	150x135x65	175x165x190	20-80	0-30
GS0102	2.0	40	60	150	150x135x100	175x165x225	20-80	0-30
GS0202	2.0	40	120	150	150x135x100	175x165x225	20-80	0-30
GS0203	3.2	40	120	150	240x135x100	265x165x230	20-80	0-30
GS0303	3.2	40	180	150	240x135x100	265x165x230	20-80	0-30
GS0304	4.5	40	180	300	300x150x100	325x175x230	20-80	0-30
GS0306	6.5	40	180	300	300x150x150	325x175x270	20-80	0-30
GS0410	10	40	240	300	300x240x150	325x265x275	20-80	0-30
GS0610	10	40	360	300	300x240x150	325x265x275	20-80	0-30
GS0615	15	40	360	450	330x300x150	355x325x275	20-80	0-30
GS0815	15	40	480	450	330x300x150	355x325x275	20-80	0-30
GS0915	15	40	540	450	330x300x150	355x325x275	20-80	0-30
GS0822	22	40	480	600	500x300x150	530x330x280	20-80	0-30
GS1022	22	40	600	600	500x300x150	530x330x280	20-80	0-30
GS1222	22	40	720	600	500x300x150	530x330x280	20-80	0-30
GS1422	22	40	840	600	500x300x150	530x330x280	20-80	0-30
GS1522	22	40	900	600	500x300x150	530x330x280	20-80	0-30
GS1030	30	40	600	600	500x300x200	530x330x330	20-80	0-30
GS1230	30	40	720	600	500x300x200	530x330x330	20-80	0-30
GS1430	30	40	840	600	500x300x200	530x330x330	20-80	0-30
GS1530	30	40	900	600	500x300x200	530x330x330	20-80	0-30

S Series Digital model (Timer & Heater & Degas) operation steps



Check the machine

Plug the machine to be sure the machine will be in a state of electricity after start the power switch on the back.

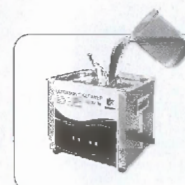

Uncover the lip and fill the tank with tap water.

No water operation will damage the machine, please put the cleaning objects in the water. Fully immerse the part to be cleaned in the water but do not over the MAX line.



Time setting

Default time is "05 00". Press "TIME+" once to increase time by 1minute. Press "TIME-" once to decrease time by 1minute. Press "ON/OFF" after time setting and ultrasonic indicator light on, machine works follow the setting parameter and stop work when time counts down to "00 00", then indicator light off.



Temperature

Temperature display panel shows "50" and actual detergent temp. Press "TEMP+" to increase setting temperature by 1°C, press "TEMP-" to decrease setting temperature by 1°C. Press "ON/OFF" to start heating function and heater indicator on, actual temperature number will blazing by 1Hz Frequency during heating process, heater indicator light will off when temperature reaches the setting temperature.



Degas

In the downtime, press "UTRASONIC" button 3 times within 3 seconds to enter the intermittent degassing mode, then the machine will work for 6 seconds and stop for 3 seconds to run alternately until the time end.



Finish work

Please do unplug the machine first after working, and then clean the tank. Be sure no liquid is in the tank if need long-term storage.



LEA ESTA INFORMACION DE SEGURIDAD

Información sobre seguridad y Recomendaciones para el uso

HOJA DE SEGURIDAD PARA RUEDAS ABRASIVAS

PRECAUCIÓN: La falta de observación de todas las precauciones le puede causar serios daños. **IMPORTANTE:** Todos los operadores deben leer y entender completamente a fondo esta información antes de usar este producto.

Advertencia: durante la operación normal de una rueda abrasiva, los materiales removidos, tales como rebabas, cascarillas, suciedad, escoria y otros residuos, serán despedidos de la rueda con fuerza considerable. Por otro lado, el paño abrasivo de la rueda se desintegra con el uso y emite partículas abrasivas y fibras.

Por esta razón existe la posibilidad de lesiones serias tanto para el operador de la rueda abrasiva como para otras personas presentes en la zona de trabajo (posiblemente hasta una distancia de 50 pies - 15 m. - o más). Con el fin de protegerse contra este peligro, los operadores y demás personas cercanas deben usar GAFAS DE SEGURIDAD o MASCARAS FACIALES COMPLETAS SOBRE GAFAS PROTECTORAS CON PANTALLAS LATERALES, junto con ROPA DE PROTECCION y EQUIPO DE RESPIRACION.

Se deben observar las instrucciones para el operador y las de seguridad además de todas las precauciones de seguridad comunes. Con ello se reducirá tanto la posibilidad como la severidad de lesiones físicas.

RESUMEN - REQUISITOS DE SEGURIDAD PARA RUEDAS ABRASIVAS

- 1. Gafas de seguridad:** TODOS LOS OPERADORES Y OTRAS PERSONAS PRESENTES EN LA ZONA de operación de ruedas abrasivas DEBEN USAR gafas de seguridad o máscaras faciales completas sobre gafas protectoras con pantallas laterales. Observe los requisitos de ANSI Z87.1, "Protección de los ojos y la cara durante el trabajo".
- 2. Cubiertas:** Mantenga las cubiertas de la máquina en su lugar. Observe la norma de seguridad ANSI B7.1.
- 3. Velocidad:** Observe las restricciones de velocidad marcadas en las ruedas abrasivas o en las cajas y etiquetas, o impresas en los manuales correspondientes. Las siglas "MSFS" significan "Maximum Safe Free Speed" o Velocidad Libre Máxima de Seguridad, es decir, girando sin aplicar trabajo (en R.P.M.). Por razones de seguridad, la "MSFS" para ruedas abrasivas montadas se especifica con el eje insertado en la herramienta de forma que la proyección del cubo de la rueda sobre el mandril o portabroca no exceda 1/4" (6 mm). No se recomienda el uso de prolongadores a menos que se reduzca la velocidad en forma correspondiente.
- 4. Equipo de protección:** Si la posibilidad de lesiones existe, se debe usar ropa y otros equipos de protección que prevengan las lesiones. Se recomienda el uso de un aparato respirador cuando la operación de la rueda abrasiva pueda producir materiales dañinos, tales como polvo de madera, partículas de metales tóxicos, etc.

Las ruedas abrasivas tienen ciertas advertencias, recomendaciones y restricciones de uso relacionadas con la seguridad marcadas directamente en el producto. No es siempre posible incluir toda la información de seguridad en la rueda. Por lo tanto, el operador DEBE LEER y OBSERVAR no sólo las instrucciones marcadas en el producto sino también las incluidas en la caja o impresas sobre ésta. El operador también debe referirse a la información de seguridad y de operación impresa en el catálogo de la rueda y en los manuales.

Condiciones de la rueda abrasiva: NO PERMITA CONDICIONES QUE NO SEAN SEGURAS. Bolilleros gastados, ejes doblados, aplicaciones poco usuales, abuso del operador o un uso no apropiado son considerados condiciones inseguras que pueden conducir a una falla de la rueda. No use o no continúe usando una rueda abrasiva que no esté funcionando correctamente (v.gr. que esté desbalanceada, etc.), esto aumenta la posibilidad de fallas y el peligro de lesiones. Toda falla o condición insegura debe ser evaluada y corregida.

⚠ ADVERTENCIA: Este producto puede exponerlo a sustancias químicas, incluyendo Dióxido de titanio, que el estado de California considera que causa cáncer. Para obtener más información, vaya a www.P65Warnings.ca.gov



LISEZ CES CONSIGNES DE SECURITE

Mesures de sécurité et conseils d'utilisation

CONSIGNES DE SECURITE POUR LES ROUES ABRASIVES

AVERTISSEMENT: De Sérieuse blessures pourrait avoir lieu à défaut de lire attentivement les précautions de sécurité **IMPORTANT :** Tout utilisateur/trices doivent lire et comprendre pleinement cette information avant l'utilisation de ce produit.

Attention: lors du fonctionnement normal des roues abrasives, les matériaux retirés, tels que les ébarbures, les écailles, la saleté, les scories de soudage ou autres résidus, sauteront avec une force considérable. Le tissu abrasif de la roue s'usera également au fil de l'utilisation, laissant échapper des grains et des fibres de tissu.

L'utilisateur de la brosse abrasive ainsi que les personnes se trouvant dans la zone de travail (dans un rayon de 15 mètres ou plus) sont exposés à des risques de blessures sérieuses. Afin de se protéger, ils doivent porter des LUNETTES DE SECURITE ou UN MASQUE FACIAL COMPLET SUR DES LUNETTES DE PROTECTION AVEC ECRANS PROTECTEURS LATERAUX, ainsi que des VETEMENTS PROTECTEURS et un EQUIPEMENT RESPIRATOIRE.

Vous devez observer toutes les instructions de fonctionnement et les consignes de sûreté ainsi que les pratiques de sécurité standard visant à réduire les risques de blessures ou leur gravité.

RÉCAPITULATIF - CONSIGNES DE SÉCURITÉ POUR L'EMPLOI DE ROUES ABRASIVES

- 1. Lunettes de sécurité:** des lunettes de sécurité ou un écran facial complet porté sur des lunettes protectrices munies d'écrans protecteurs latéraux DOIVENT ETRE PORTES par tous les UTILISATEURS de roues abrasives et toutes les PERSONNES SITUÉES DANS LA ZONE avoisinante. Agissez conformément aux spécifications ANSI Z87.1 "Occupational Eye and Face Protection" (Protection des yeux et du visage dans le cadre du travail).
- 2. Dispositifs de protection:** maintenez tous les dispositifs de protection de la machine en place. Agissez conformément aux normes de sécurité ANSI B7.1.
- 3. Vitesses:** observez toutes les limites de vitesse indiquées sur les roues abrasives, les boîtes et les étiquettes, ou imprimées dans les documents pertinents. "MFSF" signifie "vitesse maximum libre de sécurité" (tr/mn) - rotation libre sans effort. Pour des raisons de sécurité, la valeur MFSF des roues abrasives montées comprend la tige insérée dans l'outil électrique avec un surplomb maximum de 6,35mm entre le moyeu de la roue et le mandrin de l'outil. Il n'est pas conseillé d'utiliser d'arbre d'extension sauf si la vitesse de la roue abrasive est réduite en conséquence.
- 4. Equipement protecteur:** un équipement et des vêtements protecteurs appropriés doivent être utilisés dans tous les cas où ils pourraient éviter des blessures. Il est également conseillé d'utiliser un respirateur lorsque des matériaux nuisibles, tels que de la poussière de bois, des particules de métaux toxiques, etc., risquent d'être produits par l'utilisation de la roue abrasive.

Certaines mises en garde, recommandations et limitations d'utilisation relatives à la sécurité figurent directement sur les roues abrasives. Cependant, comme il n'est pas toujours possible de placer toutes les consignes de sécurité sur le produit lui-même, l'utilisateur DOIT LIRE et OBSERVER toutes les instructions fournies dans ou sur la boîte du produit ainsi que celles figurant sur le produit lui-même. L'utilisateur doit également se référer aux informations de fonctionnement et de sécurité du catalogue des roues abrasives ou de tout autre document connexe.

Etat des roues abrasives: PAREZ A TOUTES LES CONDITIONS DANGEREUSES. Tout roulement usé, arbre courbé, application inhabituelle ou utilisation abusive sont considérés comme dangereux et risquent d'entraîner une défaillance de la roue abrasive. N'utilisez pas celle-ci ou ne continuez pas de l'utiliser si elle ne fonctionne pas correctement (qu'elle n'est pas équilibrée, etc.), car ceci augmente la possibilité d'une défaillance et le risque de blessures. Toute défaillance ou condition dangereuse doit être examinée et éliminée.

⚠ AVERTISSEMENT: Ce produit peut vous exposer à des produits chimiques qui contiennent du dioxyde de titane, reconnu par l'État de la Californie comme causant le cancer. Pour de plus amples renseignements, consultez le site www.P65Warnings.ca.gov.



ULTRASONIC CLEANER

- Do not start the machine without liquid;
- With cleaning agent cleaning, better effect;
- Do not use inflammable, explosive strong corrosive cleaning agent;
- Power off after work. Dry the cleaning tank;



GRANBO

HEATER
PEAK 330 W. HI 160 W

ULTRASONIC 115 W
SEMIWAVE 52 W

8888

SET °C ACTUAL °C

0000

MINUTE SECOND

● DEGAS
● SEMIWAVE



TEMP+



TEMP-



TIME+



TIME-



HEATER



DEGAS



ULTRASONIC



SEMIWAVE

HEATER

PEAK 330 W, HI 160 W

ULTRASONIC 115 W

SEMIWAVE 52 W



GRANBO®

ULTRASONIC CLEANER

MODEL: GW0203

TIME: 0-30 min.

FREQ.: 40KHz

TANK: 3.2L

UL POWER: 60/120W

TANK SIZE: 240x135x100mm

HEATER: 150W

UNIT SIZE: 265x165x225mm

INSP. No.: QC11

INPUT: AC100V-120V/50KHz

GRANBO TECHNOLOGY INDUSTRIAL SHENZHEN CO., LTD.


WWW.GRANBOSONIC.COM

TEL:+86-755-27502790


CE FC RoHS



MADE IN CHINA

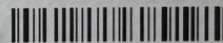

SERIAL#: 65020307690924

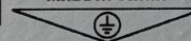

HEATER

PEAK 330 W, HI 160 W

ULTRASONIC 115 W

SEMIWAVE 52 W



GRANBO®

ULTRASONIC CLEANER

MODEL: GW0203

TIME: 0-30 min.

FREQ.: 40KHz

TANK: 3.2L

UL.POWER: 60/120W

TANK SIZE: 240x135x100mm

HEATER: 150W

UNIT SIZE: 265x165x225mm

INSP. No.: QC11

INPUT: AC100V-120V/50KHz

GRANBO TECHNOLOGY INDUSTRIAL SHENZHEN CO., LTD.

WWW.GRANBOSONIC.COM

TEL:+86-755-27502790

CE FC RoHS



MADE IN CHINA



SERIAL#: 65020307690924





DOUBLE SIDE: 45x30 CM
PREMIUM PVC: 5 LAYERS
ON
keep it flat and avoid being bent for a long time
in the service life. Mat must not be used with
liquid contact with corrosive liquids
MADE IN CHINA




11 10 9 8 7 6 5 4 3 2 1 0 1 2 3 4 5 6 7 8 9 10 11



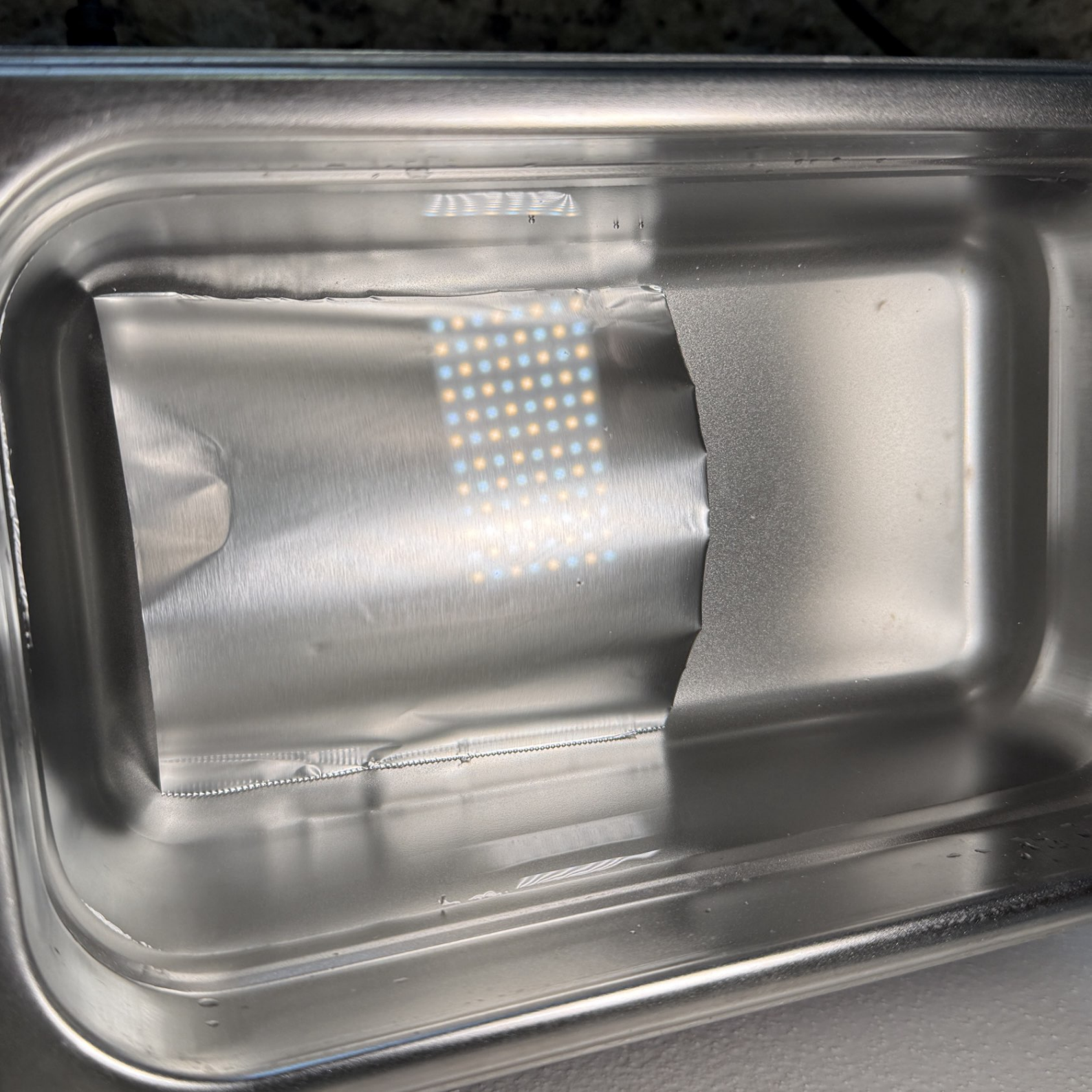

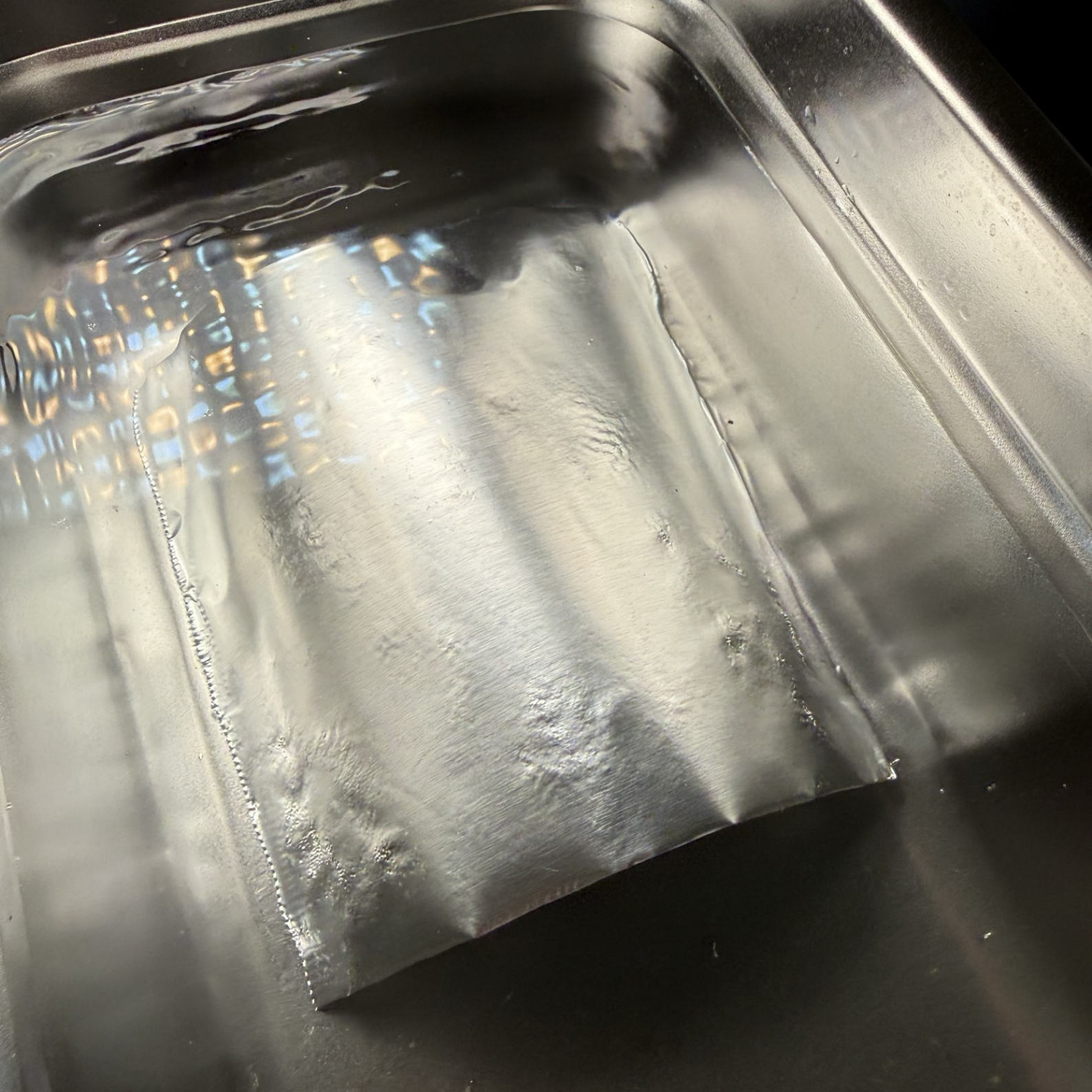


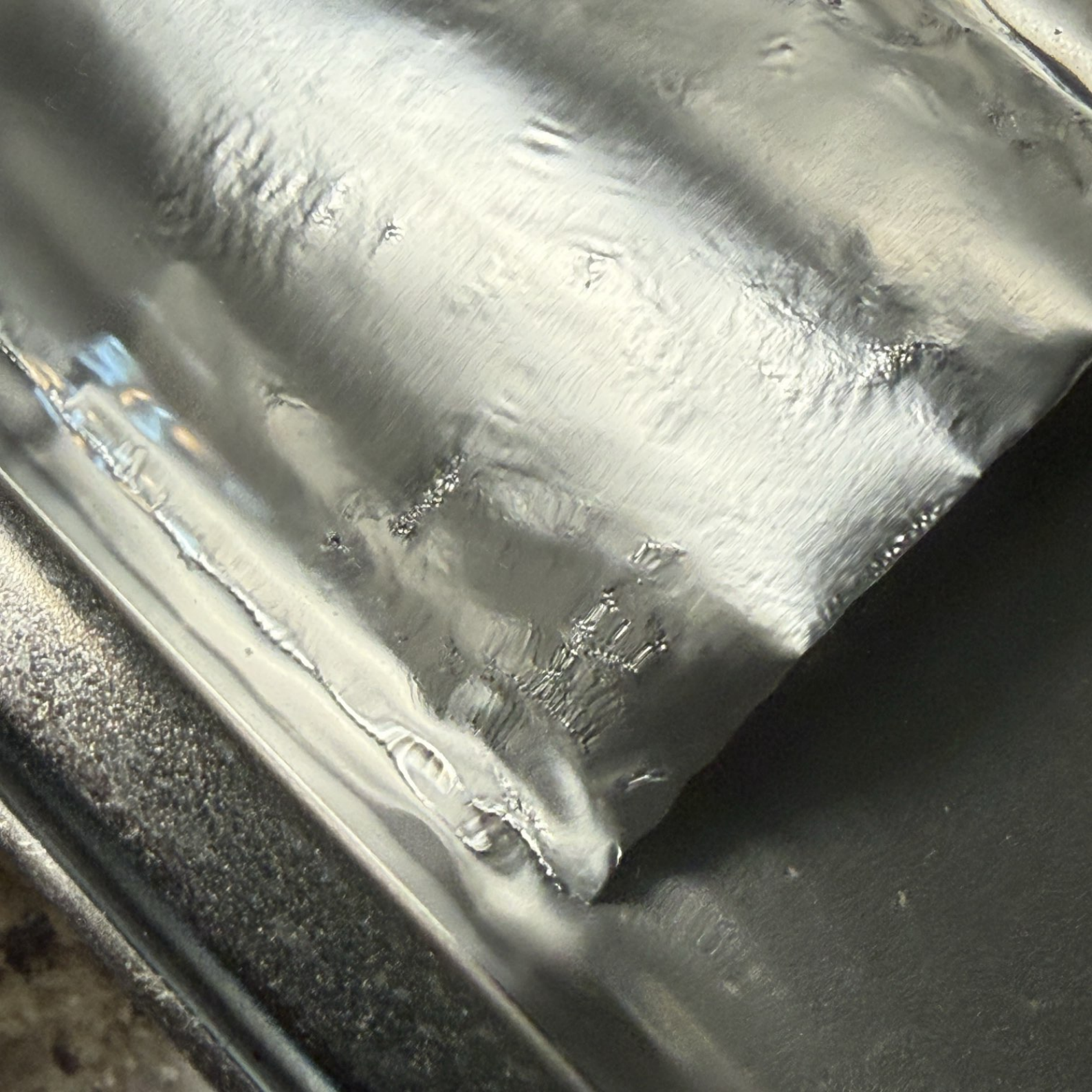



semiwave
60 sec

virgin

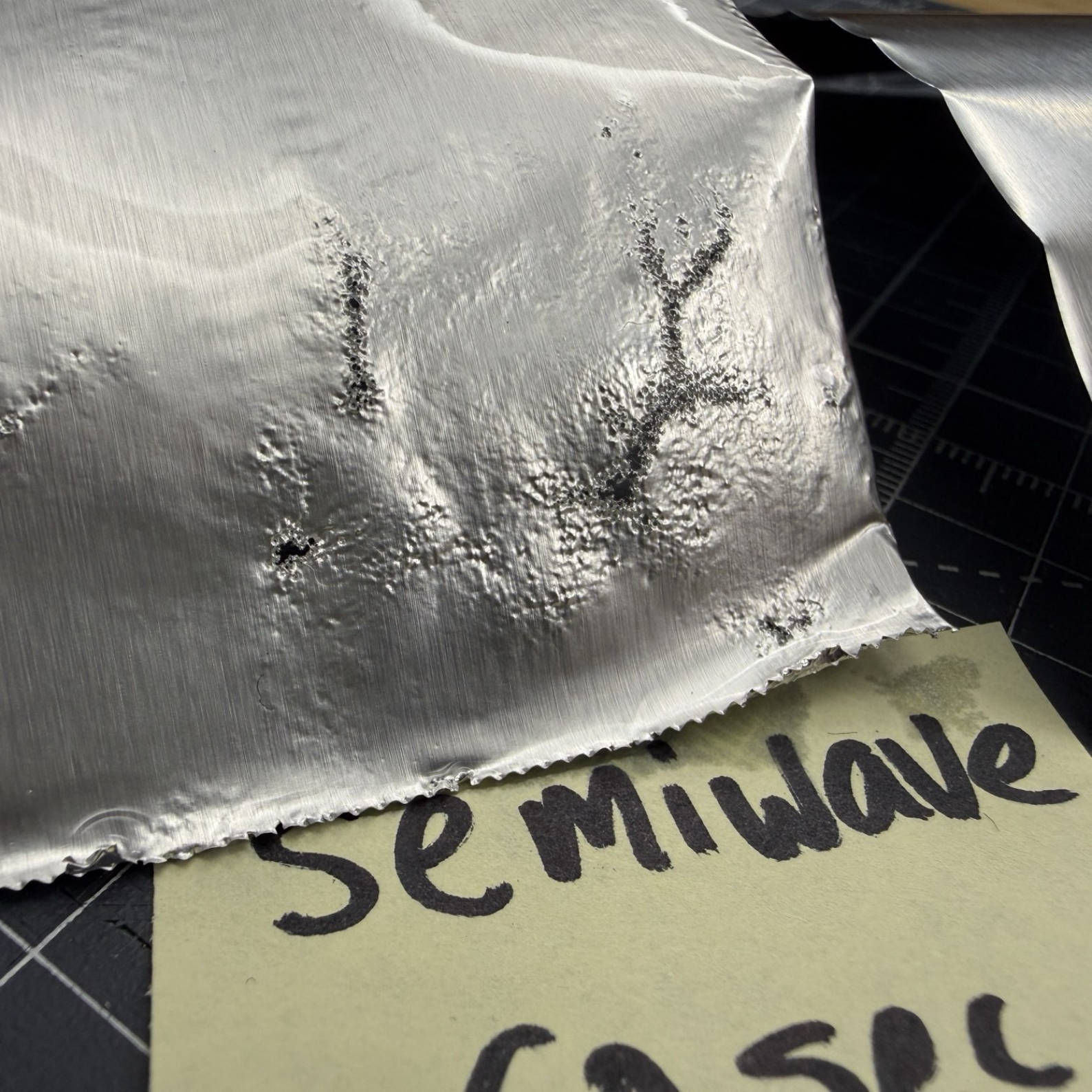

Semiwave

cas

